



# **F202 POLYMER** \* **ANTITHROMBOGENIC COATINGS**

## **REDUCTION OF CLOT FORMATION ON F202-COATED ELECTROPOLISHED STAINLESS STEEL**

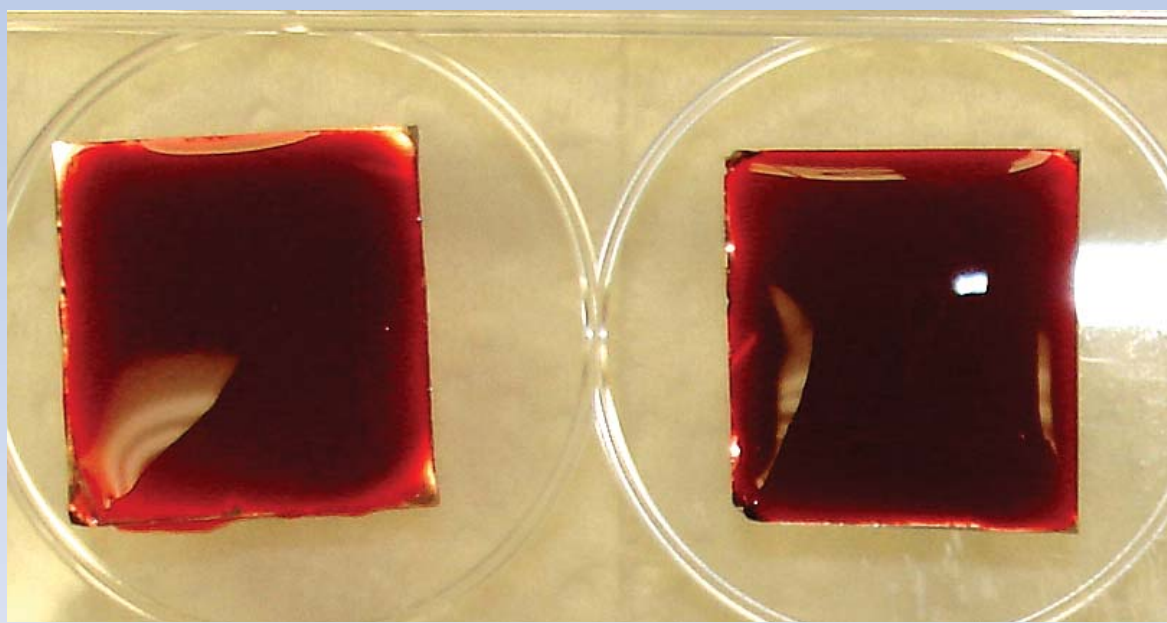


**Uncoated Control**

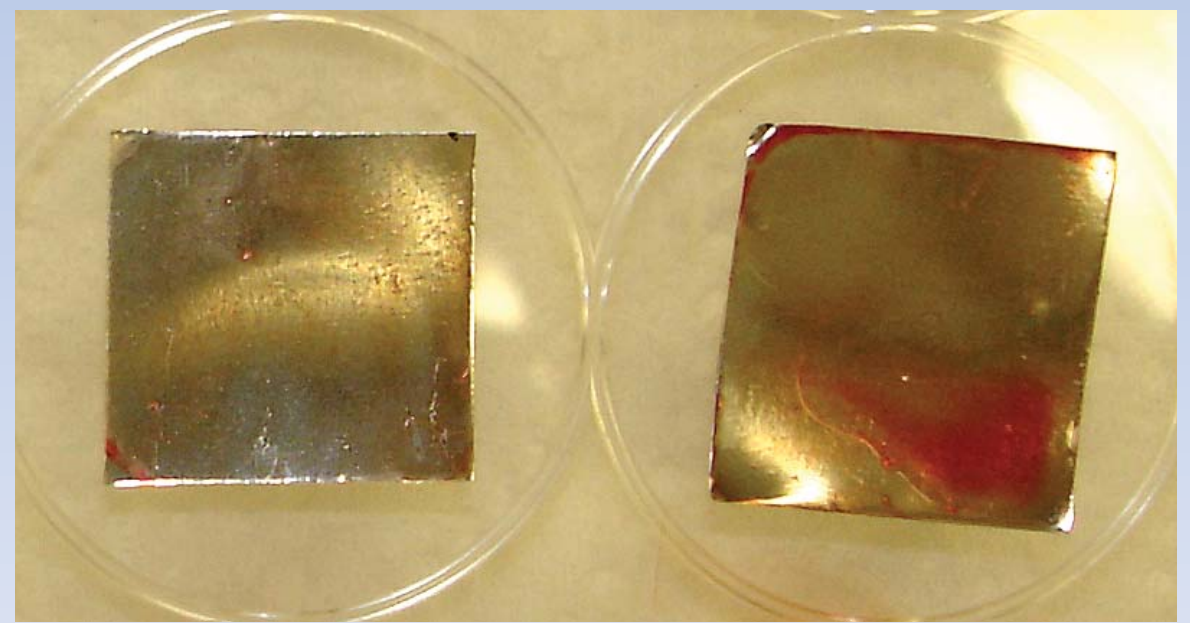


**F202-Coated**

**After extensive PBS leaching and incubations in citrated whole blood (4 x 30 minutes), recalcified whole blood was added (60 minutes), surfaces were rinsed in distilled water and examined for clot formation.**

## **ABSENCE OF CLOT FORMATION ON F200 and F202-COATED POLYURETHANE TUBING**

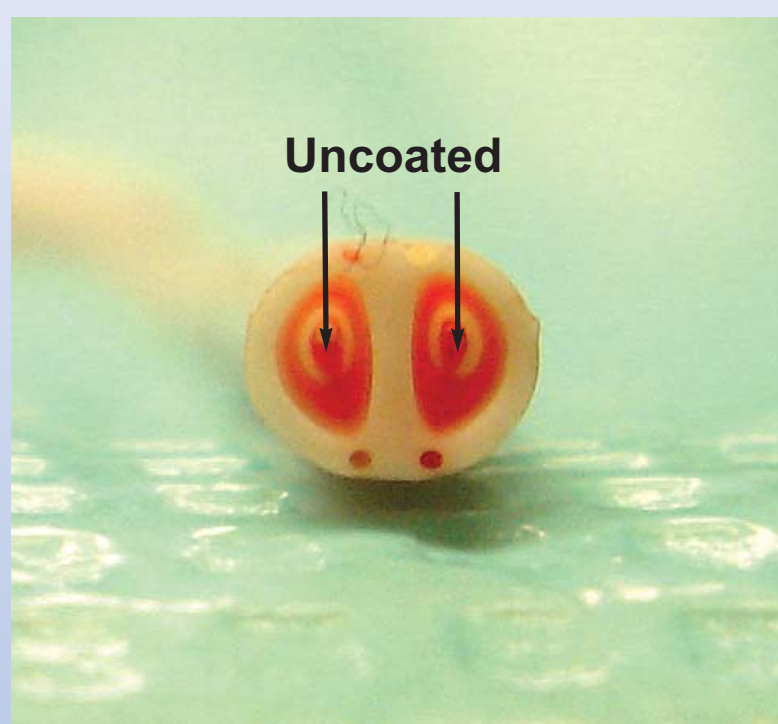

**Uncoated Control**

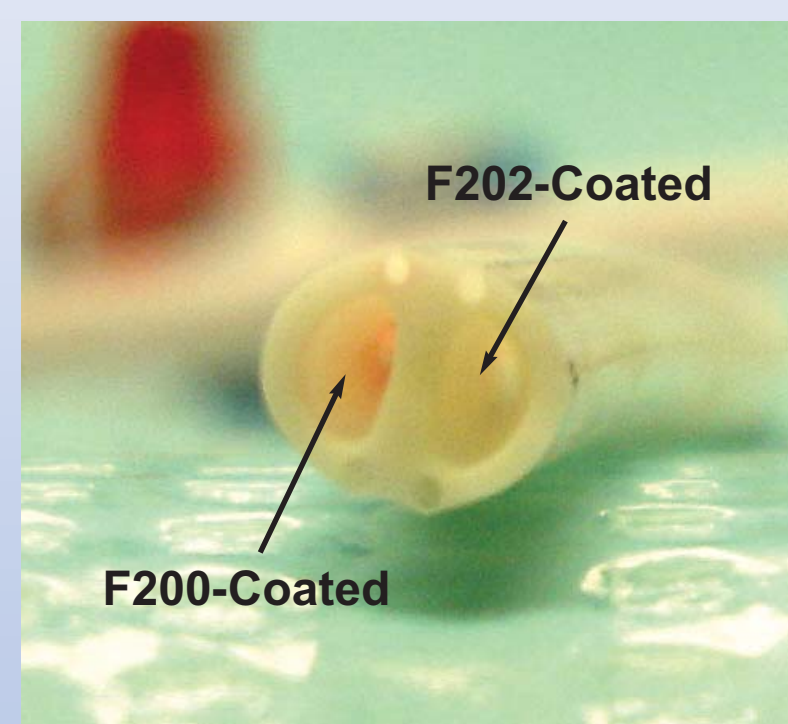

**Coated**

**\*Hydromer F202 Polymer is F200 plus bonded heparin**